

**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**

Esenyurt Firuzköy Bulvarı No:29 34325 Avcılar
İstanbul/ TÜRKİYE

TEST REPORT
DENEY RAPORU



AB-0583-T

20024487-
ing-RER

08-20

Customer name: HAYAT KİMYA SANAYİ A.Ş.
Address: SEPETLİPİNAR MAH. HAYAT CAD. NO:2 41275 YENİKÖY
BASİSKELE KOCAELİ
Buyer name: -
Contact Person: FATİH MERTTÜRK
Order No: -
Article No: GOODCARE YUMUŞAK KULAKLI MASKE
Name and identity of test item: Green non-woven mask (Claimed to be; Colour Code: Green / Green,
spunbond+meltblown nonwoven in every layer (3 layers) mask)
The date of receipt of test item: 15.07.2020
Re-submitted/re-confirmation date: -
Date of test: 15.07.2020-23.07.2020
Remarks: -
Sampling: The results given in this report belong to the received sample by vendor.
End-Use: -
Care Label: Not Specified

Number of pages of the report: 5

The Turkish Accreditation Agency (TURKAK) is signatory to the multilateral agreements of the European co-operation for the Accreditation (EA) and of the International Laboratory Accreditation (ILAC) for the Mutual recognition of test reports.

EKOTEKS LABORATUVAR ve GÖZETİM HİZMETLERİ A.Ş. accredited by TÜRKAK under registration number [AB-0583-T] for ISO 17025:2017 as test laboratory.

The test and/or measurement results, the uncertainties (if applicable) with confidence probability and test methods are given on the following pages which are part of this report.

Seal

Date
05.08.2020

Customer Representative
Hatice ACARALP

Head of Testing Laboratory
Sevim A. RAZAK

05.08.2020

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**


AB-0583-T

20024487-
ing-RER

08-20

REQUIRED TESTS	RESULT	COMMENTS
PHYSICAL PROPERTIES TESTS		
Breathability (Differential Pressure)	P	
MICROBIOLOGICAL TEST		
Bacterial Filtration Efficiency (BFE)	P	
Microbial Cleanliness (Bioburden)	P	
Blood Splash Resistance ⁽¹⁾	P	Type IIR
P: Pass F: Fail R: Refer to retailer technologist. Test results were evaluated according to EN 14683:2019+AC:2019 limit values (1) This report was reissued to add this test result.		

REMARK: Original samples are kept for 3 months and all technical records are kept for 5 years unless otherwise specified. If requested, measurement uncertainty will be reported. But unless otherwise specified, measurement uncertainty is not considered while stating compliance with specification or limit values. The reported uncertainty is based on a standard uncertainty multiplied by a coverage factor k=2, providing a level of confidence of approximately 95 %. Tests marked (*) in this report are not included in the accreditation schedule.

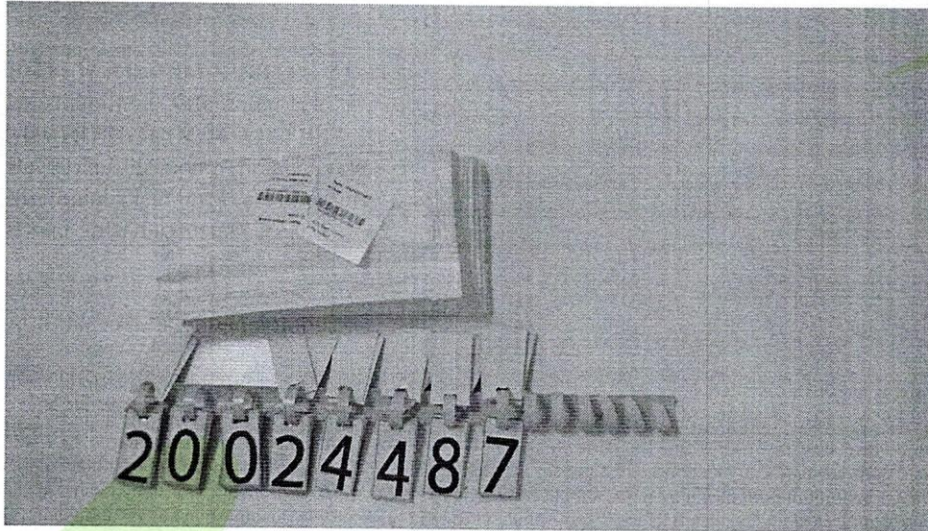

Note: The report issued on 23.07.2020 with the report number 20024487 was withdrawn and replaced with the report 20024487-RER issued on 05.08.2020 additional information of Name and identity of test item due to vendor request

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

TEST RESULT

BREATHABILITY (Differential Pressure)

Test Metodu: EN 14683:2019+AC :2019 (TS EN 14683+AC:2019) EK-C (*)

Test Condition (21 ± 5) °C ve (85 ± 5) % relative humidity, 4 hrs

Test area is 25 mm in diameter , 5 different sample was taken

Adjusted airflow is 8 l/min. The differential pressure is read directly using a differential pressure manometer .

SAMPLE	DIFFERENTIAL PRESSURE RESULT	REQUIREMENT
1	28.0 Pa/cm ²	< 40 Pa/cm ² Type I and Type II mask
2	30.8 Pa/cm ²	
3	28.8 Pa/cm ²	
4	27.0 Pa/cm ²	
5	29.8 Pa/cm ²	
Average Result	28.8 Pa/cm ²	

TEST RESULT

Medical face masks - Requirements and test methods EN 14683:2019+AC:2019 (TS EN 14683+AC:2019)

BACTERIAL FILTRATION EFFICIENCY (BFE)

Test Metodu: EN 14683:2019+AC :2019 (TS EN 14683+AC:2019) EK-B (*)

A specimen of the mask material is clamped between a impactor and an aerosol chamber. An aerosol of *Staphylococcus aureus* is introduced into the aerosol chamber and drawn through the mask material and the impactor under vacuum. The bacterial filtration efficiency of the mask is given by the number of colony forming units passing through the medical face mask material expressed as a percentage of the number of colony forming units present in the challenge aerosol.

Test Flow Rate	28,3 L/min
Total Test Flow Time	2 minute
Sample Sizes	20x20 cm ²
Test Condition	(21 ± 5) °C and (85 ± 5) % relative humidity, 4 hours
Test Microorganism	<i>Staphylococcus aureus</i> ATCC 6538
Bacterial concentration (cfu/ ml)	5x10 ⁵ cfu/ ml
incubation conditions	24 hour, 35°C ± 2°C
Positive control sample average of number of Bacteria (C)	2.4x10 ³ cfu/ ml
Mean particle size (MPS)	2.9 µm

RESULTS			Requirement BFE (%)
Number of Test Sample	Test Sample (T) Number of Bacteria (cfu)	Bacterial Filtration Efficiency (% B)	
1	30	%98.5	Type I ≥95
2	34	%98.3	Type II ≥98
3	36	%98.2	
4	34	%98.3	
5	25	%98.8	

cfu: Colony-forming unit

$$B = (C - T) / C \times 100$$

%B: Bacterial Filtration Efficiency

C: is the mean of the total plate counts for the two positive control runs

T: is the total plate count for the test specimen

TEST RESULT

MICROBIAL CLEANLINESS (Bioburden)

Test Metod: EN 14683:2019+AC :2019 (TS EN 14683+AC:2019) EK-D (*)
EN ISO 11737-1:2018 /TS EN ISO 11737-1 :2018 (*)

5 sample were taken. The sample is weighted and put in extraction liquid after shaking well (250 rpm, 5 min), inoculated on the suitable agar.
The plates are incubated for 3 days at 30 ± 1 °C for 72 hours, and 7 days at (20 to 25) °C for TSA and SDA plates respectively. Total microorganisms counts are calculated.

	RESULTS	REQUIREMENT
Microbial cleanliness (cfu/g)	10 cfu/g	≤ 30 cfu/g Type I and Type II mask

*cfu= Colony forming unit.

BLOOD SPLASH RESİSTANCE

Test Metod: EN 14683:2019+AC :2019 (Clause 5.2.4) the resistance of the medical face mask to penetration

ISO 22609 :2004 Clothing for protection against infectious agents — Medical face masks — Test method for resistance against penetration by synthetic blood (fixed volume, horizontally projected)

Test Condition (21 ± 5) °C ve (85 ± 5) % relative humidity, 4 hrs
5 different sample was taken

	<u>SPLASH RESISTANCE PRESSURE (kPa)</u>	<u>RESULTS</u>	<u>REQUIREMENT</u>
1	>21.3 kPa	PASS	≥ 16 kPa Type IIR mask
2	>21.3 kPa	PASS	
3	>21.3 kPa	PASS	
4	>21.3 kPa	PASS	
5	>21.3 kPa	PASS	
Average Result	>21.3 kPa	PASS	